


Bipolar Stimulation Catheter

TYPE: SOFT TIP

The specific design of the soft and flexible tip prevents damage of vessel walls. Thus, it can be recommended particularly for patients with delicate vessel structures.


For temporary use in case of emergency, for pre- or postoperative stimulations, for interim support during pacemaker exchange and for atrial and ventricular system monitoring.

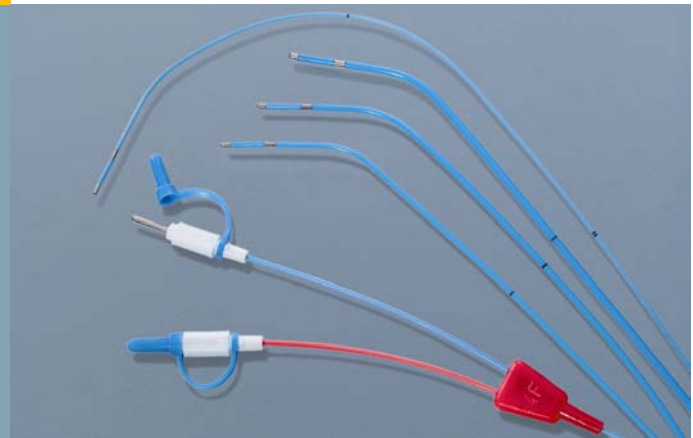

Temporary Stimulation Catheter with soft tip

FEATURES:

- unique design: particular soft, flexible tip with solid catheter body
- excellent torque control
- pliable material
- increased safety by shrouded contact pin
- indication at distal contact pin
- radiopaque
- length markings

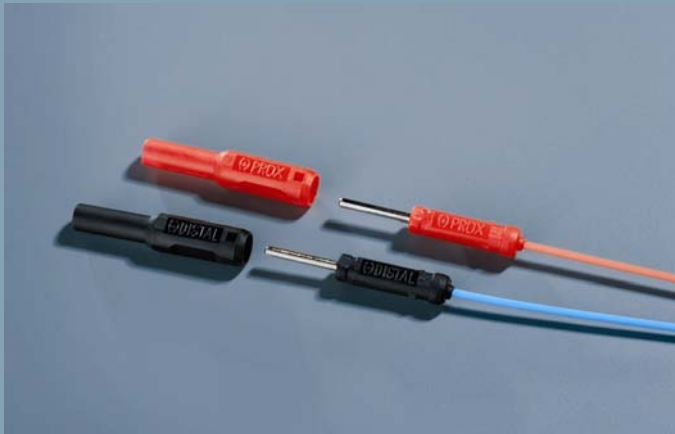

Shrouded pin adapter

Description	Ref.No.	Size	Electrode Surface	Introducer / FEP Cannula
Temporary bipolar stimulation catheter with soft tip Contents: Stimulation catheter with FEP Cannula, standard pin and shrouded pin adapter	25.TPL5	5 F / 1.7 mm	15 mm ²	5 F
	25.TPL6	6 F / 2.0 mm	19 mm ²	6 F
	25.TPL70	7 F / 2.3 mm	22 mm ²	

STERILE EO

SPECIFICATIONS:

- Material:
 - Catheter Tube:..... PUR, radiopaque
 - Electrodes:..... Stainless steel
- Dimensions:
 - Length of flexible tip:..... 70 mm
 - Usable length:..... 110 cm
 - Number of electrodes:..... 2
 - Electrode tip:..... hemispherical
 - Electrode distance:..... 10 mm
- Packaging:.....double packed, sterile
- Packing Unit:..... 5 EA / Box