


Saugermed®

Neurosurgical Brain Pads

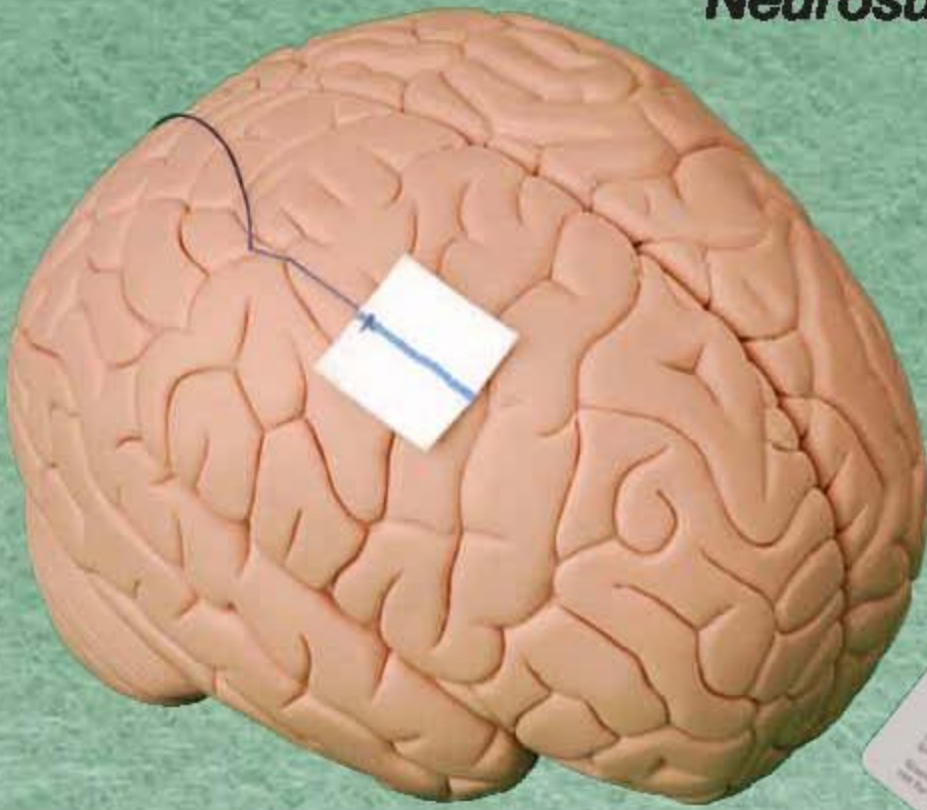

Neurosurgical Products

ISO 13485:2003

CE 1984




SAUGERMED®

SURGICAL BRAIN PADS

Within the context of Class III for neurosurgical operations Special produced haemostatic featured pads

Various options for different use of purpose 9 different featured and structured pad produced from 9 various raw material

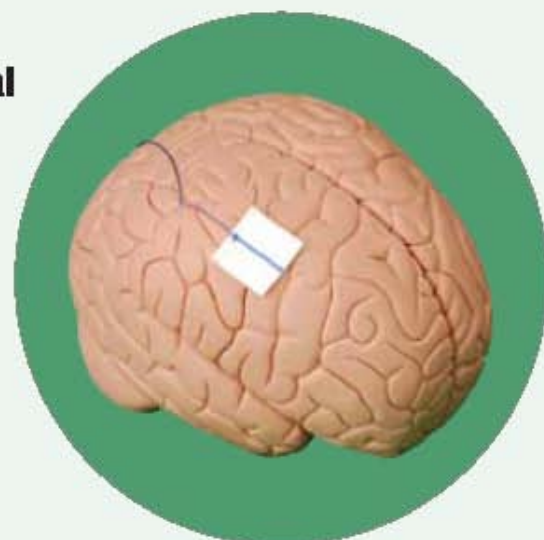

Saugerméd *thinpad*

- Absorbation capacity 8 times of its weight
- Very thin, smooth and silky structured pad



Saugerméd *plus*

- Absorbation capacity 8 times of its weight
- Thin, closely structured pad



Saugerméd *thincot*

- Absorbation capacity 9 times of its weight
- Thin, smooth structured, laminated surface pad



Saugerméd *laminat*

- Absorbation capacity 10 times of its weight
- Standard thickness, smooth structured, laminated surface pad



Saugerméd *soft*

- Absorbation capacity 10 times of its weight
- Standard thickness, closely structured, laminated surface pad



Saugerméd *ultra soft*

- Absorbation capacity 15 times of its weight
- Smooth, delicate structured pad



Saugerméd *cotton soft*

- Absorbation capacity 16 times of its weight
- Very smooth, delicate structured pad



Saugerméd *purecotton*

- Absorbation capacity 17 times of its weight
- Very thin, very light structured cotton pad



Saugerméd *purecotton soft*

- Absorbation capacity 23 times of its weight
- Very soft, very absorbative, closely structured cotton pad

* Absorbation capacity of the Neurosurgical Brain Pads certified by the accredited laboratory report approved by the Notify Body

COMMON FEATURES

Atraumatic

Produced from hydrofiber structured organic and inorganic materials in nonwoven form
Atraumatic and soft.

Unfuzzy and unfeathered, does not leave loofah

Enables safe use

All the pads consists of radio-pack blue strap which is X-ray detectable

Pls. see figure 1

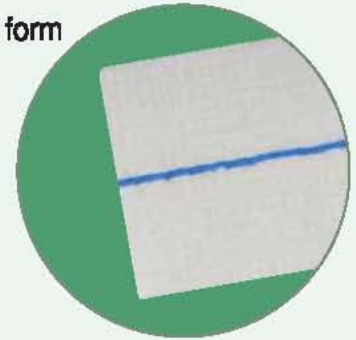


Figure 1

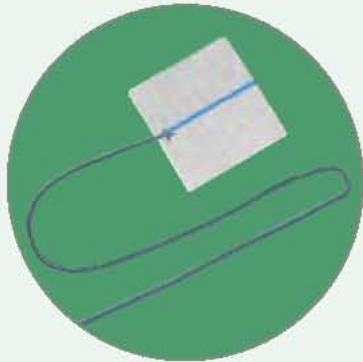


Figure 2

Also, all the pads have x-ray detectable atraumatic featured

Trace rope in case of smooth touch

By its trace rope it can be easy to identify the place of the pad and take out

Pls. see figure 2

***Patent rights of the x-ray detectable trace rope of Saugermed protected by the certification of the Turkish Patent Institute as Utility Model**

To avoid any distention Trace rope looped with a special technics

Pls. see figure 3

***Patent rights of the special loop protected by the certification of the Turkish Patent Institute as Utility Model**

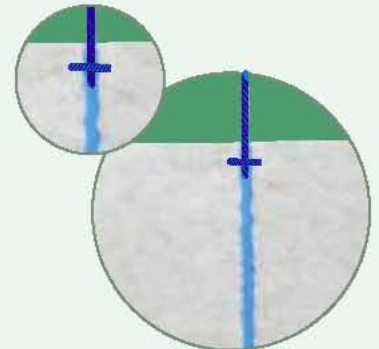


Figure 3



Figure 4

To recognize the number of used pads during the operation it has a Special Counter Card 1 to 10

Pls. see figure 4

To prevent the pads from deformation it consists of a protective belt

Pls. see figure 5



Figure 5

Sterile and dirtiness protected


Gamma sterilized

Tvyek made of sterile pouch prevents unbounded tear and can only be open from its opening part

It has an anticontamination bag in the sterile pouch.

There are 16 different dimensions of each pad made of 9 various raw material

Pls.see Order Table

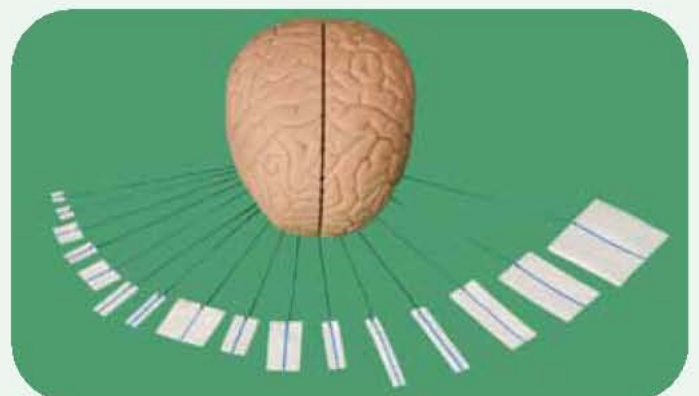

ORDER TABLE

Dimension (mm)	Piece Pouchette	Pouchette Box	THINPAD	PLUS	THINCOT	LAMINATE	SOFT	ULTRA SOFT	COTTON SOFT	PURE COTTON	PURE COTTON SOFT
10 x 10	10/1	24/1	BPTP 1010	BPPLS 1010	BPTC 1010	BPLMN 1010	BPSFT 1010	BPUSFT 1010	BPCSFT 1010	BPPC 1010	BPPCS 1010
10 x 40	10/1	24/1	BPTP 1040	BPPLS 1040	BPTC 1040	BPLMN 1040	BPSFT 1040	BPUSFT 1040	BPCSFT 1040	BPPC 1040	BPPCS 1040
10 x 75	10/1	24/1	BPTP 1075	BPPLS 1075	BPTC 1075	BPLMN 1075	BPSFT 1075	BPUSFT 1075	BPCSFT 1075	BPPC 1075	BPPCS 1075
12.5 x 12.5	10/1	24/1	BPTP 1212	BPPLS 1212	BPTC 1212	BPLMN 1212	BPSFT 1212	BPUSFT 1212	BPCSFT 1212	BPPC 1212	BPPCS 1212
12.5 x 25	10/1	24/1	BPTP 1225	BPPLS 1225	BPTC 1225	BPLMN 1225	BPSFT 1225	BPUSFT 1225	BPCSFT 1225	BPPC 1225	BPPCS 1225
12.5 x 40	10/1	24/1	BPTP 1240	BPPLS 1240	BPTC 1240	BPLMN 1240	BPSFT 1240	BPUSFT 1240	BPCSFT 1240	BPPC 1240	BPPCS 1240
12.5 x 50	10/1	24/1	BPTP 1250	BPPLS 1250	BPTC 1250	BPLMN 1250	BPSFT 1250	BPUSFT 1250	BPCSFT 1250	BPPC 1250	BPPCS 1250
12.5 x 75	10/1	24/1	BPTP 1275	BPPLS 1275	BPTC 1275	BPLMN 1275	BPSFT 1275	BPUSFT 1275	BPCSFT 1275	BPPC 1275	BPPCS 1275
20 x 20	10/1	24/1	BPTP 2020	BPPLS 2020	BPTC 2020	BPLMN 2020	BPSFT 2020	BPUSFT 2020	BPCSFT 2020	BPPC 2020	BPPCS 2020
20 x 40	10/1	24/1	BPTP 2040	BPPLS 2040	BPTC 2040	BPLMN 2040	BPSFT 2040	BPUSFT 2040	BPCSFT 2040	BPPC 2040	BPPCS 2040
25 x 25	10/1	24/1	BPTP 2525	BPPLS 2525	BPTC 2525	BPLMN 2525	BPSFT 2525	BPUSFT 2525	BPCSFT 2525	BPPC 2525	BPPCS 2525
25 x 50	10/1	24/1	BPTP 2550	BPPLS 2550	BPTC 2550	BPLMN 2550	BPSFT 2550	BPUSFT 2550	BPCSFT 2550	BPPC 2550	BPPCS 2550
25 x 75	10/1	24/1	BPTP 2575	BPPLS 2575	BPTC 2575	BPLMN 2575	BPSFT 2575	BPUSFT 2575	BPCSFT 2575	BPPC 2575	BPPCS 2575
40 x 40	10/1	24/1	BPTP 4040	BPPLS 4040	BPTC 4040	BPLMN 4040	BPSFT 4040	BPUSFT 4040	BPCSFT 4040	BPPC 4040	BPPCS 4040
40 x 70	10/1	24/1	BPTP 4070	BPPLS 4070	BPTC 4070	BPLMN 4070	BPSFT 4070	BPUSFT 4070	BPCSFT 4070	BPPC 4070	BPPCS 4070
70 x 70	10/1	24/1	BPTP 7070	BPPLS 7070	BPTC 7070	BPLMN 7070	BPSFT 7070	BPUSFT 7070	BPCSFT 7070	BPPC 7070	BPPCS 7070

1/1 SCALE PRODUCT TABLE (mm)

